

# Bulk, Sculpt & Cure with Tetric EvoCeram® Bulk Fill

**NEW  
PRODUCT**

No additional viscosities.  
No additional layers.  
No additional equipment.  
*No Compromise!*



**SPECIAL OFFER!**

Valid November 1, 2011 – March 31, 2012

**BUY: 2 Refills GET: 1 Refill FREE!\***



Introducing  
**Tetric EvoCeram® Bulk Fill**  
NANO-HYBRID COMPOSITE RESTORATIVE

100% CUSTOMER SATISFACTION  
GUARANTEED!

[ivoclarvivadent.com](http://ivoclarvivadent.com)

Call us toll free at 1-800-533-6825 in the U.S.  
© 2011 Ivoclar Vivadent, Inc. Tetric EvoCeram is a registered trademark of Ivoclar Vivadent.

**ivoclar  
vivadent**  
passion vision innovation

# Introducing Tetric EvoCeram® Bulk Fill

a unique nano-hybrid composite material developed specifically for the fast, efficient, "bulk placement" of direct posterior composite restorations. **Simply Bulk, Sculpt and Cure.**



## BULK

The material's smooth consistency facilitates excellent adaptation to cavity walls eliminating the need for flowable liners or additional dispensing equipment.



## SCULPT

The well-balanced filler composition results in excellent physical properties allowing the material to be contoured immediately after placement without the need for a final "capping" layer.




## CURE

The patented "Polymerization Booster" provides a fast, complete cure.


## Polymerization Booster

Tetric EvoCeram Bulk Fill features an innovative photo-initiator system offering greater reactivity to curing lights as compared to traditional initiators, allowing 4mm increments of composite material to be completely cured in just 5 seconds (eg, bluephase 20i @ 2,000 mW/cm<sup>2</sup>).

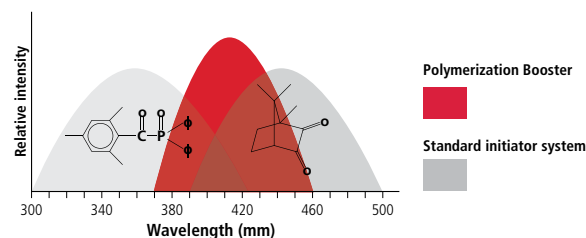


## Light Sensitivity Inhibitor

A patented additive incorporated into the photo-initiator system delays the effects of ambient light (e.g. operating light) maintaining a sculptable consistency until polymerization is desired, providing longer working time as compared to other bulk fill materials.

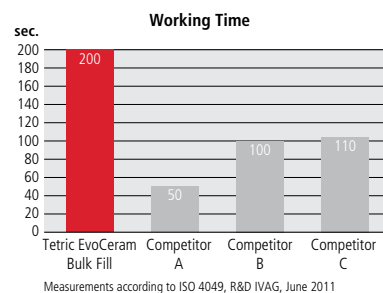

## Shrinkage Stress Reliever

Patented Shrinkage Stress Relievers integrated into the filler composition provide the lowest combination of volumetric shrinkage and shrinkage stress for superior marginal integrity.

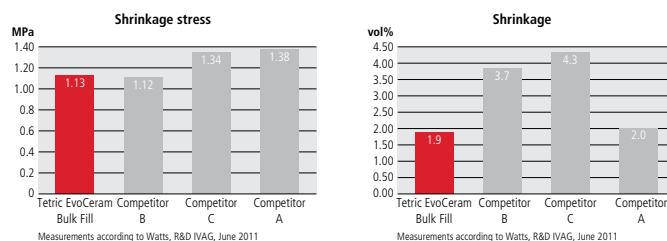

\*Not a trademark of Ivoclar Vivadent, Inc.

## 3 Universal Shades

The three universal shades; IVA, IVB and IVW (white shade for fast deciduous fillings or light-colored teeth,) demonstrate enamel-like translucency allowing restorations to blend seamlessly with natural dentition.



\*Offer valid: 11/1/11 – 3/31/12. Qualifying purchases must appear on a single invoice. Ivoclar Vivadent reserves the right to modify or cancel this promotion at any time. This special cannot be combined with any other offer or rebate program.

**FREE GOODS REDEMPTION:** Mail or fax a copy of your Patterson invoice no later than 4/30/12 to Ivoclar Vivadent, Inc., Attn: Clinical Promotions, 175 Pineview Drive, Amherst NY 14228 – FAX (716) 691-2285. Indicate PROMO CODE: 612 and shade preference of free product(s). Shade IVA will be sent if no preference is indicated. Cannot mix & match delivery forms.

100% CUSTOMER SATISFACTION  
GUARANTEED!

[ivoclarvivadent.com](http://ivoclarvivadent.com)

Call us toll free at 1-800-533-6825 in the U.S.

© 2011 Ivoclar Vivadent, Inc. Tetric EvoCeram is a registered trademark of Ivoclar Vivadent.

ivoclar  
vivadent®  
passion vision innovation