

DENTAL LIGHT-CURE BONDING AGENT

CLEARFIL™ S³ BOND PLUS

[CLEARFIL TRI-S BOND PLUS]

Stronger than Before, Faster than Ever and Simple as Always.

CLEARFIL S³ BOND PLUS is a light-cure, self-etching, one-step, fluoride-releasing bonding agent. Based on the same progressive adhesive technologies from the gold standard; **CLEARFIL SE BOND** you can expect minimal postoperative sensitivity.

With its one-bottle liquid, **CLEARFIL S³ BOND PLUS** dramatically simplifies bonding procedures. Along with Kuraray's **MDP** technology, a new catalyst has been adopted; thereby, providing higher bond strength with improved working time.

CLEARFIL S³ BOND PLUS has minimal technique sensitivity, and is indicated for light-cure composites, root surface treatment and cavity sealing. **CLEARFIL S³ BOND PLUS** is exceptional for use in pedodontic restorative procedures where speed and comfort are extremely critical.

CLEARFIL S³ BOND PLUS when utilized with our complementary component **CLEARFIL DC CORE PLUS** can achieve dual-cure core build-up/post cementation without an additional bottle.



kuraray

DENTAL LIGHT-CURE BONDING AGENT

CLEARFIL™ S³ BOND PLUS

[CLEARFIL TRI-S BOND PLUS]

One-Step Self-Etch: No-Shaking, No-Mixing, No-Scrubbing, No Multiple Applications

SAME FEATURES AS CLEARFIL S³ BOND

- Swift and Simple: No-Shaking, No-Mixing, No Scrubbing, No Multiple Applications
- **MDP Technology** Provides a Water-Resistant Chemical Bond to Enamel and Dentin without Etching
- Perfect for Root Surface Treatment and Cavity Sealing
- Virtually No Post-Operative Sensitivity

“PLUS” FEATURES

- Fluoride-Releasing
- New Integrated Catalyst for Advanced Bond Strength
- Post Cementation and Dual-cure Core Build-Ups using **CLEARFIL DC CORE PLUS**

KURARAY'S MDP TECHNOLOGY:



CLEARFIL S³ BOND PLUS contains the same original monomer **MDP** created by Kuraray, which is the key for reliable adhesion to tooth structure and restorations.

MDP Technology enables **CLEARFIL S³ BOND PLUS** to chemically interlock with the calcium in both enamel and dentin, creating a highly durable water-proof seal.

MDP has been utilized in the well-recognized and established dental adhesive and resin cement:

- **CLEARFIL SE BOND**
- **PANAVIA**

QUICK & SIMPLE APPLICATION:

1. Apply bond for 10 sec.
2. Dry with mild pressure air flow for 5 sec.*
3. Light-cure for 10 sec. (curing time by Halogen/LED)

* Bond contains ethanol, which could lower the bond strength if it is not dried completely.




ORDERING INFORMATION:

ITEM#	PRODUCT NAME	CONTENTS
#2890KA	CLEARFIL S ³ BOND PLUS Set	1 Bond (4 ml) 1 Applicator brush (fine 50 pcs) 1 Dispensing dish 1 Light-blocking plate
#2893KA	CLEARFIL S ³ BOND PLUS VALUE Pack	3 Bond (4 ml)
#2895KA	CLEARFIL S ³ BOND PLUS Single Dose Pack	50 Unit-dose 0.1 ml/tip
#2896KA	CLEARFIL S ³ BOND PLUS Single Dose VALUE Pack	100 Unit-dose 0.1 ml/tip

Comparison of CLEARFIL S³ BOND PLUS and other Leading Bonding Agents

Tensile Bond Strength to:
Enamel and Dentin after
4000 Thermal-Cycles

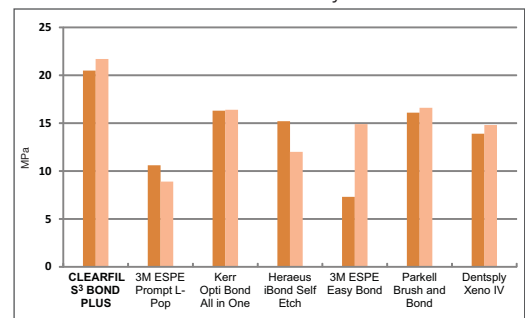

■ Enamel ■ Dentin *Source: KURARAY MEDICAL INC.

RELATED PRODUCT:



CLEARFIL DC CORE PLUS

is a dual-cure, radiopaque core build-up material supplied in a convenient automix delivery system, which is perfect for post cementation. With **CLEARFIL S³ BOND PLUS** it's fast and easy; providing exceptional bond strength to dentin.

Warranty: KURARAY MEDICAL INC. will replace any product that is proved to be defective. Neither KURARAY MEDICAL INC. nor the distributor accepts liability for any loss or damage, direct, consequential or special, arising out of the application or use of or the inability to use these products. Before using, the user should determine the suitability of the products for the intended use and the user assumes all risk and liability whatsoever in connection therein. CLEARFIL™ S³ BOND PLUS and KURARAY MEDICAL INC. products are available only through Kuraray Authorized dealers.

CLEARFIL™ and PANAVIA™ are trademarks of KURARAY Co., Ltd.